


# NEUATION

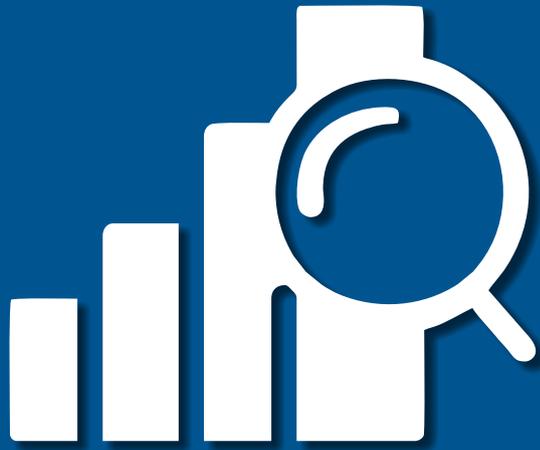

VARIABLE UTILITY  
PRECISE RESULTS

iTherm D150-1  
(Dry Bath)

## iTherm D150-1 (Dry Bath)

When using a digital thermal bath the problems are numerous. From uneven heating to fine tuning the device to suit your requirements, more often than not using it feels like a terrible task. Low temperature range, corrosive surface, and non-uniform heating are some of the most common issues faced by users which create more problems for them rather than simplifying the whole process.

Backed by strong R&D into customer expectations and their feedback, we at Neuation came up with the iTherm D150-1 which combines exceptional performance with a maintenance-free long life. Designed to provide a superior user-experience, the device provides uniform heating across the blocks to ensure reliable performance. Additionally, numerous user-friendly features like variable parameter setting, auto temperature safety cut etc. ensure that a technician's work in the lab is made easier without expending much effort or time.



## Product Features



Dry Heating Block With Temperature Range From Ambient +5°C To 120°C.



$\leq 0.5^{\circ}\text{C}$  Of Accuracy And  $\pm 0.5^{\circ}\text{C}$  Of Uniformity Across The Blocks



Lid For Protection And Heat Preservation



Rapid Heating And Uniform Temperature Across The Block Throughout Its Operation.

Composite Raw Material Body For Higher Chemical Resistance

**NO  
RUST**

Variety Of Heating Block For Specific Tubes



Can Be Used As Water Bath Or Bead Bath



Large And Clear Display Shows Critical Parameters And All Set Values In One Glance Without Any Toggle.




## Product Features




Wide Tube Capacity Range  
With Convenient Lifting Gear  
For Cleaning And Disinfecting.

NTC1000K probe for  
accurate measurement  
of temperature

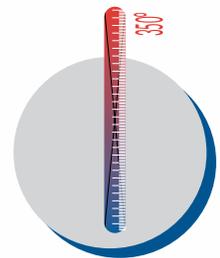

## HEATING BLOCKS

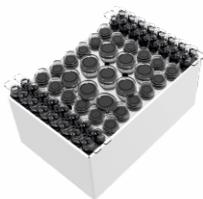

**Hybrid Block**  
(9 x 1.5/2 ml, 20 x 0.5 ml  
32 x 0.2/0.1 ml PCR tube)



**50 ml Tube Block**  
(5 x 50ml)



**15 ml Tube Block**  
(12 x 15 ml)



**1.5/2 ml Block**  
(28 x 1.5/2 ml)



**Bead Block**  
(Liquid: 300ml (water)  
Beads: 1.5kg)






**Heating Block 28x2 ml**  
(10x75mm  
Round Bottom Tube)

## Ordering information

| Specifications         | iTherm D150-1            |
|------------------------|--------------------------|
| Temperature Range      | Ambient +5 to 120°C      |
| Temperature Uniformity | ± 0.5°C                  |
| Temperature Stability  | ± 0.5°C                  |
| Set-up Block Material  | Aluminium                |
| Number Of Blocks       | 1                        |
| Display                | Digital                  |
| Timer                  | 19Hrs 59 Mins & Infinite |
| Weight                 | 1.7Kg                    |
| Dimension              | 190x87x82                |
| Protection Class       | IP 21                    |
| Input Voltage          | 100-240 VAC 50/60Hz      |
| Output Voltage         | DC 24 V 4A               |
| Power Consumption      | 85W                      |

# Ordering information

| iTherm D150-1 |                       |                      |
|---------------|-----------------------|----------------------|
| Product Name  | Cat #                 | Details              |
| iTherm D150-1 | DB02-1111 EU/US/UK/AU | Single Block Drybath |
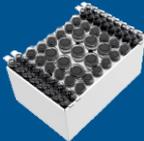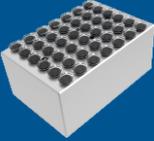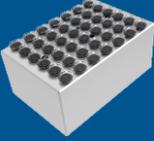
| Accessories                                    |  |         |   |
|--|--|---------|---|
| Product Name                                   | Capacity                                       | Cat #   | Product   |
| Hybrid Block                                   | 9x1.5/2ml<br>20x0.5ml<br>32x0.2/0.1ml PCR Tube | DOB1-B1 |    |
| 1.5 ml MCT Block                               | 28 x 1.5 ml                                    | DOB1-B2 |    |
| 2 ml MCT Block                                 | 28 x 2 ml                                      | DB02-B6 |   |
| 15ml Tube Block                                | 12x15ml  | DB02-B3 |  |
| 50ml Tube Block                                | 5x50ml   | DB02-B4 |  |
| Heating Block<br>28x2ml<br>(Round Bottom Tube) | 2ml Tube<br>(10x75mm)                          | DB02-B5 |  |
| Bead Block                                     | Liquid: 300ml (Water)<br>Beads: 1.5kg          | DB02-BB |  |
| Beads  | DB02-B   | DB02-B  | -   |

# NEUATION

www.neuation.com

AN ISO 9001 : 2015 & ISO 13485 : 2016 CERTIFIED MFG. COMPANY

\*We also execute OEM projects. Specifications can be changed without notice for quality improvement.